


QA3™

Patient Stretcher

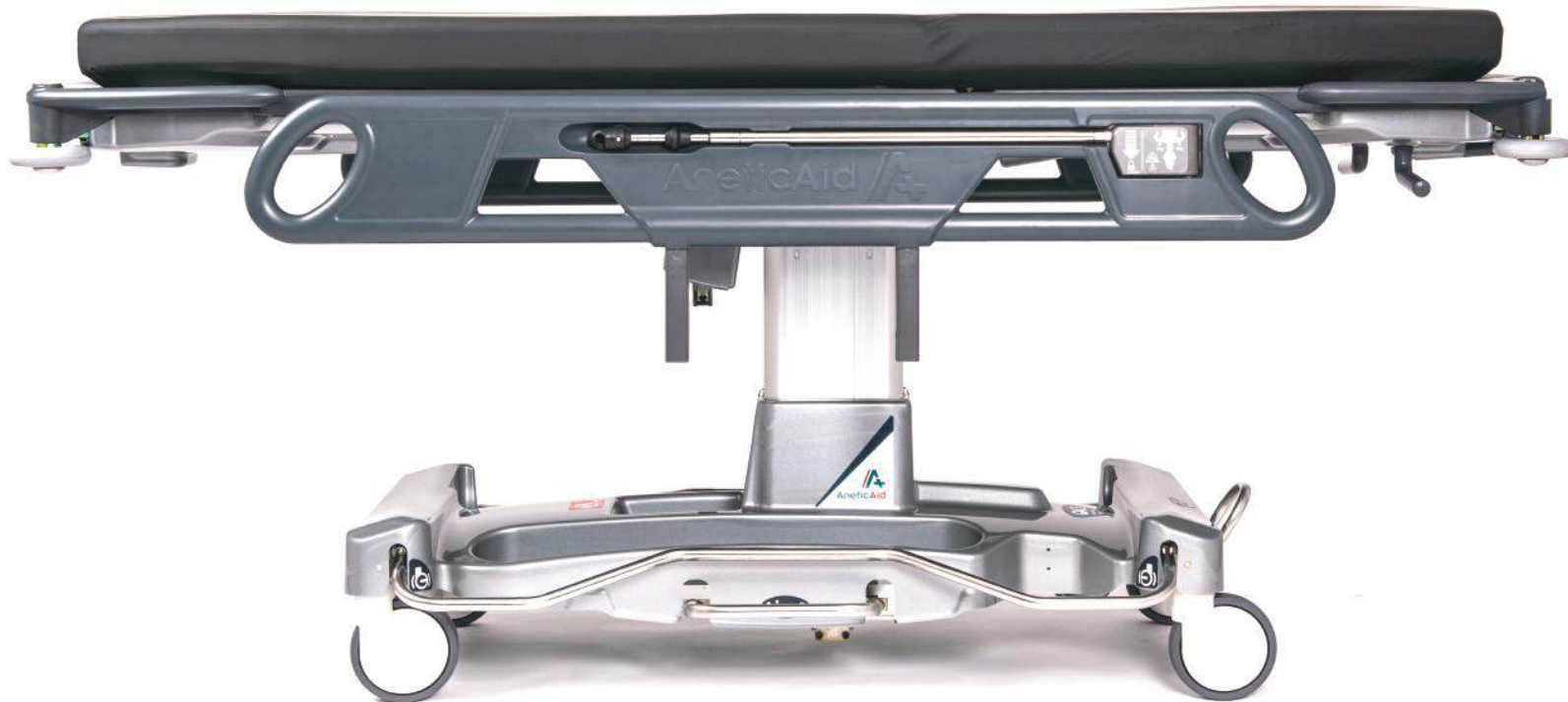

The complete solution for **patient transport**

AneticAid® 

QA3 PATIENT STRETCHER SYSTEM

QA3 **|||||DRIVE**™ Patient Stretcher

**Our flagship stretcher -
with inbuilt motorised drive.**

Reduce your manual handling liabilities further with the compact, light weight and supremely manoeuvrable QA3™ DRIVE - the complete solution for patient transport.

Ideal for patient transport, examination and intubation in theatre and clinical areas, the QA3™ DRIVE is robust and has maximum patient comfort built in. The stretcher comes with exceptional low height capabilities and new pressure care mattress - now 5% more resistant to chlorine and other disinfectants.

Featuring inbuilt motorised drive, operated by innovative touch-sensitive pushing handles, the QA3™ Drive can be effortlessly controlled and travel up to gradients of 6°.

Featuring:

- Inbuilt motorised drive
- Innovative drive handles
- 5th wheel steering
- Long-life battery



Description	Code
QA3™ DRIVE Patient Stretcher	21112



Three models to choose from



QA3™ Treatment Stretcher

Our brand new powered model

Supreme comfort, optimal positioning and optimal access for a wide range of clinical treatments - the QA3™ Treatment Stretcher offers both the user and the patient and improved clinical journey.

Description	Code
QA3™ Treatment Stretcher	21116



QA3™ Patient Stretcher

Standard features - our flagship model

Ideal for patient transport, examination and intubation in theatres and clinical areas. Supremely manoeuvrable, compact, lightweight and robust with maximum patient comfort built in.

Description	Code
QA3™ Patient Stretcher	21110



DESIGN & INNOVATION

Combining industry knowledge with
innovative technology



MADE IN THE UK

Designed and produced at our UK manufacturing
facility using the latest engineering techniques



GLOBAL DISTRIBUTION

Manufacturer, designer and supplier of high quality
medical equipment around the world



Subject to annual servicing

Unique manufacturer-backed Lifetime Warranty Scheme

STANDARD FEATURES

Supremely manoeuvrable, compact, lightweight and robust with maximum patient comfort built in. The QA3 Patient Stretcher is ideal for patient transport, examination and intubation in theatre and any other clinical area.

With the backrest raised, it is little more than 1.8m long, descends lower and is narrow and shorter than nearly all other patient stretchers on the market.

Gas assisted backlift

Unique design, precision controlled, easy lift system for enhanced user and patient comfort.



Exceptional patient comfort

Low height of 475mm (excluding mattress) & store away side rails allow for easier mount/dismount, particularly for frail and elderly patients. Also including a versatile positioning headrest, wide patient platform and pressure care mattress for comfort.

Trendelenburg tilt

Extremely smooth and controlled mechanism, even at lowest height.

Pushing handles

Four ergonomically positioned pushing handles for supreme manoeuvrability

5th wheel steering

Our unique mechanism ensures maximum manoeuvrability and means the QA3 is easily guided by one support worker.

Simultaneous braking low rolling resistance castors

Polyurethane tyres offer low rolling resistance for high manoeuvrability: 150mm diameter central total lock castors fitted as standard, helping to offer a rigid and reliable work platform.

Vertically operated safety side rails

Store away under the mattress improving the ease of mounting and dismounting the stretcher for the patient.

Infection control

Smooth moulded surfaces and acrylic capped vac formings enclose trolley mechanics to minimise dirt traps. Anti-bacterial properties and innovative design for easy cleaning.

Moving and handling

Zero transfers mean risk is minimised for patient and practitioner alike: from the ward, through anaesthetics, to the surgical procedures, and finally to first and second recovery.

Fixed transfusion pole


Innovative quick release safety locking mechanism, integrated into the right-hand foot end side rail for easy access and storage.

Cylinder trough

Side loading for ease of access, suitable for wide range of cylinder options including F, E or CD

Hydraulic lift mechanism

Exceptional low height capability, adjusted by independently operating pedals on both sides of stretcher

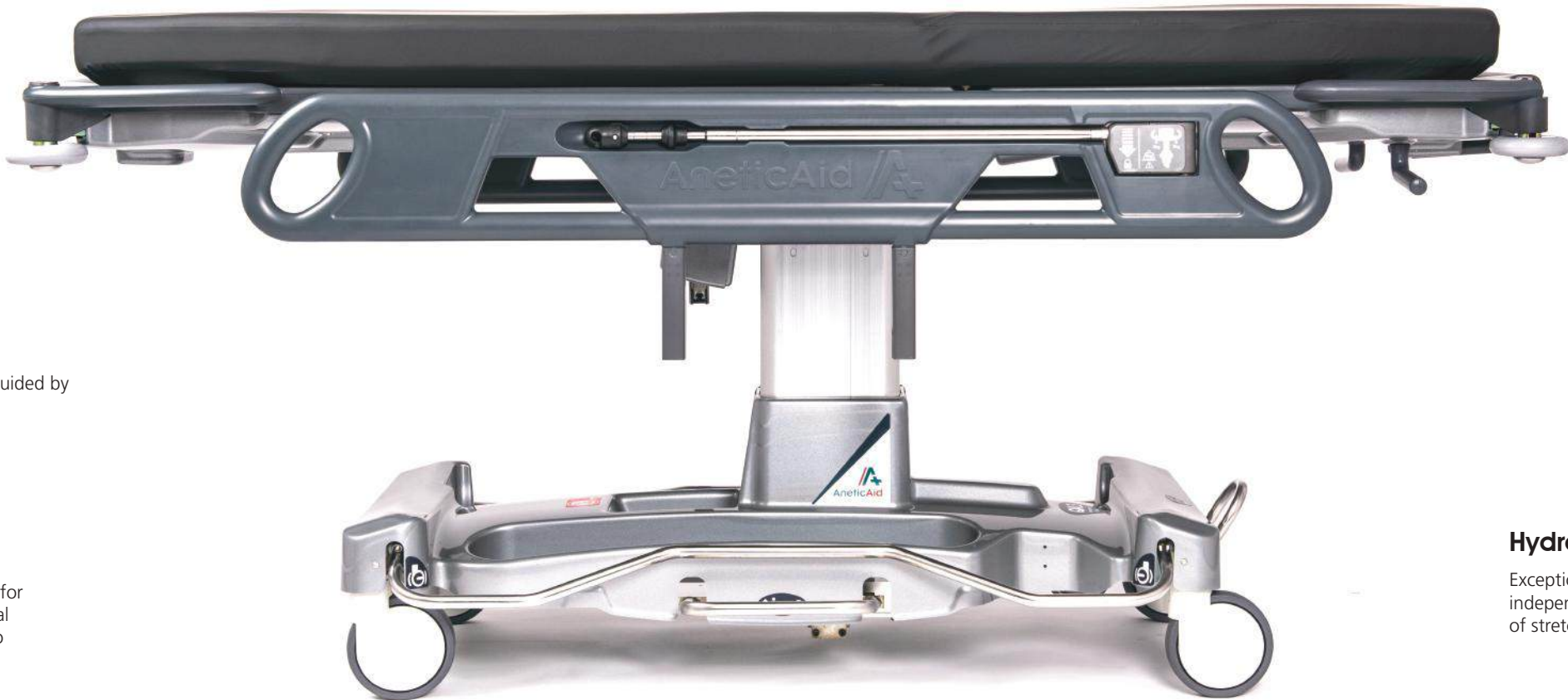

Pressure care mattress

Manufactured from temperature-sensitive gel-feel foam that moulds to the patient's body shape, relieving pressure by distributing the load as evenly as possible. Outer material is anti-microbial, water resistant, vapour permeable and latex free. Now 5% resistant to chlorine and other disinfectant compounds.



SPECIFICATIONS		Code: 21110/21112	
Minimum height	*475mm	Mattress length	2025mm
Maximum height	*775mm	Mattress width	700mm
*excludes mattress		Mattress depth	90mm
Minimum length	1700mm	Weight	128kg
Maximum length	2150mm	Backrest	0 to +90°
Overall width	805mm	Trendelenburg	0 to -12°
Castor diameter	150mm	Reverse Trendelenburg	*0 to -10°
Maximum patient weight			320kg

*NOTE - Reverse Trendelenburg enabled through a build option

Robust with maximum patient comfort built in



DESIGN & INNOVATION



MADE IN THE UK



GLOBAL DISTRIBUTION

///DRIVE™ ASSIST FEATURE

Inbuilt motorised drive offers effortless movement, combined with ultimate control.

Because people and therefore patients are getting larger, it's crucial that we take moving and handling issues into account. We want to give practitioners a stretcher on which patients can be moved and positioned as easily as possible.

The introduction of our touch-activated, motorised DRIVE feature offers a major enhancement in a stretcher that already has ease of use, practicality and patient comfort built in. The QA3™ DRIVE can also be used in traditional fifth wheel steer mode.



Effortless travel

The QA3™ DRIVE Patient Stretcher can travel up to gradients of 6°.



Inbuilt motorised drive

No need to push or pull the stretcher: our inbuilt motorised DRIVE feature linked to the stretcher's fifth wheel offers superb manoeuvrability and control



Long-life battery

Capable of 20 kilometers on a single charge, the DRIVE feature is powered by an advanced nickel metal hydride rechargeable battery.



Accurate battery readings

Highly accurate battery charge indicator is combined with a rapid charge function.



How does ///DRIVE™ ASSIST work?

QA3 ///DRIVE
Patient Stretcher

Here's the breakdown...

Firstly, how does DRIVE assist work functionally? Designed into the backrest frame of QA3 DRIVE are strain gauges that are connected to push handles.

As you push and pull the QA3 DRIVE, you flex the strain gauge frame imperceptibly, and this force is translated through its software into a torque requirement for the DRIVE assist wheel.

The greater the force needed, the greater the DRIVE assist delivered and done so in a natural, intuitive manner. You walk and a QA3 DRIVE will move at your pace — supreme control, supreme confidence.



How does it benefit you?

Secondly, how does DRIVE assist work practically?

According to health and safety executive guidelines, when operating a stretcher or patient trolley, the pushing or pulling force **should not equal or exceed 20 kilos to start a load in motion**. The pushing or pulling force should also not equal or exceed 10 kilos to keep a load in motion.



Cornering with a load exerts greater lateral forces on the body, affecting the shoulders, hips and lower back.

Which is why Anetic Aid have introduced the innovative and intuitive QA3 DRIVE Stretchers - to vastly reduce the start/stop and lateral forces whilst manoeuvring a patient transport stretcher.



DESIGN & INNOVATION



MADE IN THE UK



GLOBAL DISTRIBUTION

Motorised drive: effortless travel, ultimate control

⚡ POWERED FEATURES

Minimise risks for practitioners dealing with increasing patient weights with powered positioning: the whole stretcher platform or the backrest can be raised or lowered electronically from the panel at the foot end.

With six powered functions – platform raise and lower, back rest, KneeFlex, Profile, Trendelenburg and Reverse Trendelenburg - as well as CPR return – the QA3™ Treatment Stretcher is a truly powerful and versatile platform for transport, treatment and recovery.

Battery powered
Long-life rechargeable lithium-ion battery (NMC) with highly accurate charge indicator



Electrically Powered Functions
Easy to operate using the hand control positive feel push button controls.

SPECIFICATIONS		Code: 21116	
Minimum height	*495mm	Mattress length	2025mm
Maximum height	*875mm	Mattress width	700mm
*excludes mattress		Mattress depth	90mm
Minimum length	1700mm	Maximum patient weight	320kg
Maximum length	2150mm	Backrest	0 to +70°
Overall width	805mm	KneeFlex	0 to +30°
Castor diameter	150mm	Trendelenburg	0 to -12°
Weight	134kg	Reverse Trendelenburg	0 to -10°



Platform height



Offers smooth movement to position the stretcher at any point in its height range – 495mm to 875mm



Backrest

Positions the backrest anywhere between the vertical and horizontal for maximum patient comfort



KneeFlex

Easy elevation of the knees to offer improved comfort in a supine position



Profiling

Easy elevation of the back rest and knees at the same time to offer improved comfort during recovery



CPR Return

If a patient needs CPR treatment, there is no time to lose: the powered CPR feature returns the stretcher swiftly to a horizontal, optimal height position for CPR treatment



Foot-end Control

Even greater flexibility when repositioning the QA3™ using the foot-end panel



Trendelenburg / Reverse Trendelenburg

Extremely smooth and controlled by a powered actuator, even at lowest height



Effortless positioning at the touch of a button



DESIGN & INNOVATION



MADE IN THE UK



GLOBAL DISTRIBUTION

MORE MODELS TO CHOOSE FROM

QA3™ Powered ⚡ Emergency Stretcher

The complete solution for any emergency department - now with powered functions.

Reduce your manual handling liabilities further with the compact, light weight and supremely manoeuvrable QA3™ Powered Emergency Stretcher - the complete solution for any emergency department.

Specially designed for Emergency Departments with four powered functions and a radiotranslucent platform, the QA3™ Powered Emergency Stretcher offers exceptional cassette access, is robust and has maximum patient comfort built in. It also includes exceptional low height capability and a pressure care mattress now 5% resistant to chlorine and other disinfectants compared to the previous 1%.

Powered functions include: platform height, backrest, Trendelenburg/Reverse Trendelenburg and CPR return.

- Featuring:
- Four powered functions
 - X-ray translucent platform
 - Fixed transfusion pole
 - Exceptional patient comfort



Description	Code
QA3™ Powered Emergency Stretcher	21126



Discover the QA3 Emergency range



QA3 /// DRIVE™ Emergency Stretcher

Standard features with inbuilt DRIVE function

Featuring inbuilt motorised drive, operated by innovative touch-sensitive pushing handles, the QA3™ DRIVE Emergency Stretcher offers controlled, effortless travel up to gradients of 6°.

Description	Code
QA3™ DRIVE Emergency Stretcher	21122



QA3™ Emergency Stretcher

Standard features - our complete solution for emergency departments

Specially designed for Emergency Departments, the QA3™ Emergency Stretcher is perfect for pre-treatment and post-treatment patient transport.

Description	Code
QA3™ Emergency Stretcher	21120



QA3 /// DRIVE™ <800> Bariatric Stretcher


Our specially designed bariatric stretcher with DRIVE function

The QA3™ <800> DRIVE Bariatric Stretcher has all the benefits of the QA3™ DRIVE range, now with an extended 800mm platform - specially designed for bariatric patients.

Description	Code
QA3™ <800> DRIVE Bariatric Stretcher	21115

ACCESSORIES

Accessorise your QA3 Patient Stretcher at the time of ordering to better suit your department needs or add greater versatility to your QA3 with this range of convenient accessories.

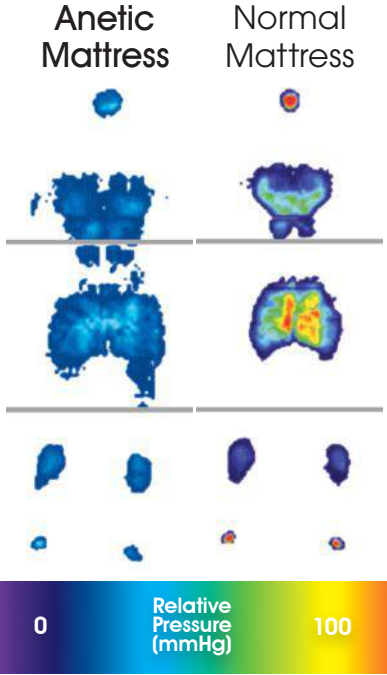

Pressure Care Mattress

Manufactured from temperature-sensitive gel-feel foam that moulds to the patient's body shape, relieving pressure by distributing the load as evenly as possible. Outer material is anti-microbial, water resistant, vapour permeable and latex free. Now 5% resistant to chlorine and other disinfectant compounds compared to the previous 1%.

Description	Code
Replacement Mattress with Waterproof Zip	21186
Replacement Mattress Cover with Waterproof Zip	21186-B1

Side-Rail Covers

Foam padded covers, available in both a grey and child-friendly print, are designed to completely encapsulate the QA3™ side rail.




Description	Code
Side Rail Cover - PADDED - GREY Material	21172
Side Rail Cover - PADDED - CHILD PRINT Material	21173

Theatre Pillow

Tired of the number of pillows heading to the laundry, reducing their availability? This superbly comfortable theatre pillow is covered in a water-impermeable material and reduces the turnaround time for pillows between uses. Simply clean & disinfect by hand between uses.



Description	Code
Theatre Pillow	10545

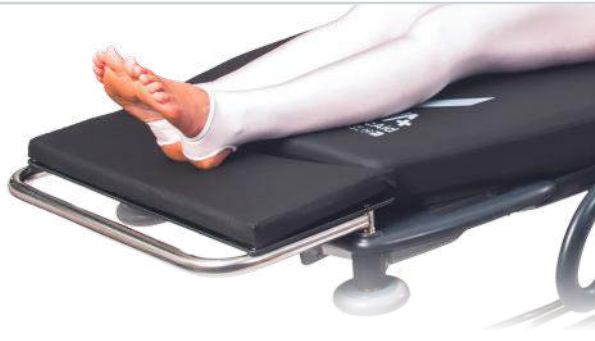

Foot-end Extension

Need extra leg room? Specifically designed for use with the QA3™, this foot-end mounted accessory adds an extra 300mm in length to the stretcher for taller patients.

Description	Code
Foot-end Extension	21153



Endoscopy Demi-frame

Peace of mind protection. This demi-frame side rail is specifically designed with Endoscopy units in mind.

Description	Code
Side Rail - Endoscopy Height Extension Demi-frame	21165



Refreshment Tray

Specially designed for use with a QA3™, this innovative shelf has retractable legs which can be deployed and safely seated onto the top of the side rails - becoming a convenient patient refreshment tray.

Description	Code
Refreshment Tray Only	21152-A1



QA3 Accessory Bracket

A new side rail accessory bracket which can be fitted to our QA3™ Stretchers, allowing for other QA3 accessories to be fitted with ease.

Description	Code
QA3 Accessory Bracket	21156

Thinking about hand surgery? Think Anetic.

We have three specialist arm and hand surgery systems fully compatible with our QA3 Accessory Bracket - all for different positions and procedures.

Try out:

- The Gardiner Hand Surgery System
- The Magna Hand Surgery Table
- The Woodlands Arm & Hand Surgery Table

If you're thinking about solutions for arm and hand surgery, then we've got you covered.



Accessorise your QA3™ to suit your department needs



BUILD OPTIONS

Customise your QA3 Patient Stretcher at the time of ordering to better suit your department needs or add greater versatility to your QA3 with this range of build options.



Reverse Trendelenburg

Enable 12 degrees of Reverse Trendelenburg to your QA3 Stretchers with this handy build option. Extremely smooth and controlled mechanism allows this to be achieved even when the stretcher is at its lowest height.



Patient's Possession Storage Container

In need of more storage space on your QA3™? Add a conveniently positioned, foot-end base mounted storage container, mounted onto retractable runners for when the container isn't needed.



Department ID Labelling

Struggling to identify a QA3™ as belonging to your department? Add your department name to your QA3 to aid equipment identification.

PLEASE NOTE: A MAXIMUM OF 22 CHARACTERS INCLUDING SPACES ARE ALLOWED.

Additional Transfusion Pole - Right Hand Side Rail Mounted



Add a second fixed transfusion / infusion stand to the right hand side rail on your QA3™ Patient Stretcher. Popular with Endoscopy departments.



Oxygen Cylinder Bracket - For CD Sized Cylinder

Add a compact and conveniently positioned bracket, retaining the position of CD sized Oxygen Cylinder brackets within the base mounted cylinder trough.



Bracketing for Emergency Suction Systems

Did you know? We can offer a compact and conveniently positioned bracket for your emergency / SOS suction system to be mounted on your QA3 Stretcher.

This build option allows your QA3 Stretcher to be suitable for a wide range of suction devices - perfect for critical procedures that require a patient's airways to be clear.

***NOTE: Oxygen & Suction Apparatus not sold by Anetic Aid.**



Subject to annual servicing




We offer the unique and popular 'Lifetime Warranty Scheme'

Take out a **AnetiCare** service package for regular and predictable maintenance costs that includes a **10 YEAR WARRANTY**.

Covers engineer's time and most parts. Contact us to learn more about what an **AnetiCare** service package can offer you.

Available via your Anetic Aid partner – find them at aneticaid.com/distributors.

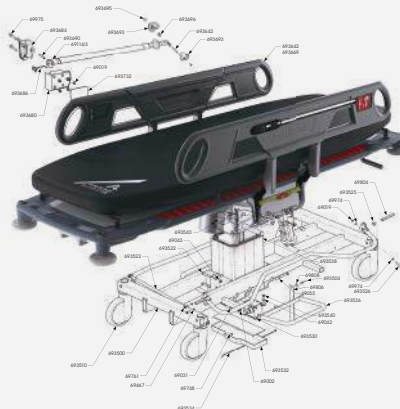

DESIGN & INNOVATION



MADE IN THE UK



GLOBAL DISTRIBUTION

Customise your QA3 with our range of build options



QA3 Ophthalmic Stretcher:

Outstanding surgical access - exceptional patient comfort.



PLUS
our other product areas



Scan this code to see our **QA4's** journey through day surgery!



SEE
ALSO

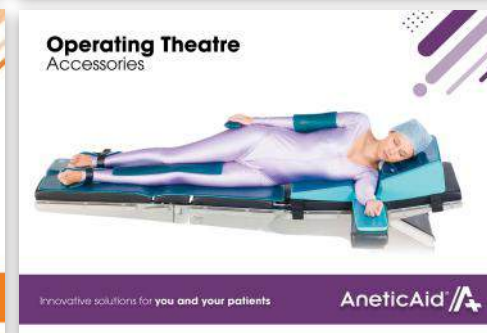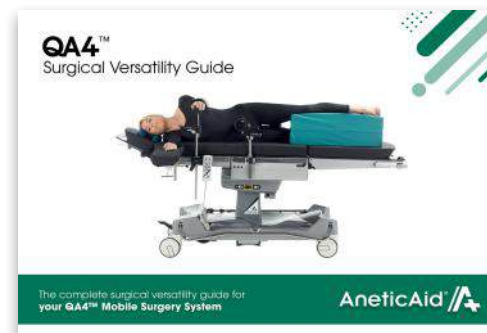

For more information on any of our products or service contracts, call:
sales@aneticaid.com +44 (0)1943 878647 aneticaid.com



DESIGN & INNOVATION



MADE IN THE UK



GLOBAL DISTRIBUTION

Innovative **Medical** Technology - **Practically Applied**



Anetic Aid Ltd, Sapper Jordan Rossi Park, Otley Road, Baildon BD17 7AX, UK