

MC Series

**LAMIDEY
NOURY**
MEDICAL



Your reliable partner
for electrosurgery



Monopolar / Bipolar currents

A real haemostatic cut with a white coagulation without charring and resection in liquid without bleeding.



Endosection

Haemostatic cutting current dedicated to polypectomy



Saline cut bipolar

Dedicated to hystero-resection procedure



Micro-coelio bipolar

Specific precise coagulation without charring nor sticking dedicated to neurosurgery, spinal and laparoscopic surgery

UNITS RANGE

MC 2 / 200 W




MC 3 / 370 W



MC 4 / 370 W

2 Monopolar outlets



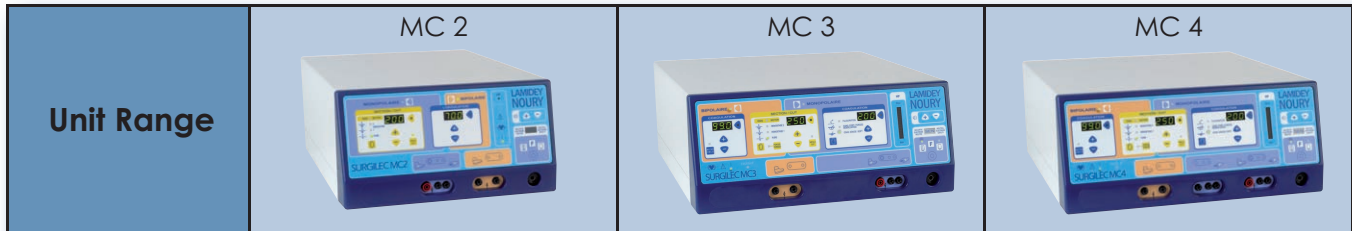

TECHNICAL SPECIFICATIONS

Main supply	230 V 50 - 60 Hz / 500 VA	230 V 50 - 60 Hz /1200 VA	230 V 50 - 60 Hz /1200 VA
Size and weight	34.5x19x37 cm / 9 kg	42x19x43cm / 11 Kg	42x19x34cm/ 11Kg
HF nominal frequency	550 KHz	450 KHz	450 KHz
Frequency modulation	20 KHz	21 KHz	21 KHz

OUTPUT POWER

Maximum power	200 W	400 W	400 W
Monopolar cut	200 W	400 W	400 W
Monopolar coag	80 W	200 W	200 W
Bipolar cut	—	250 W	250 W
Bipolar coag	70 W	100 W	100 W

ERGONOMY

Power display	Digital display 3 digits	
Outlets	1 monopolar / 1 bipolar	2 monopolars / 1 bipolar

SAFETIES

S.E.N Neutral Electrode Safety
CPA Automatic power control
Automatic limitation of the RF voltage
EXCESS HF power safety
High Frequency leakage limitation
Defibrillation shock safety
Short circuit safety
Class I, type CF



Lamidey Noury Medical

Z.A. Les Godets 3 rue des Petits Ruisseaux
F 91370 Verrières le Buisson France
Tél. : +33 (0)1 69 20 69 69
Fax : +33 (0)1 60 13 97 47
info@lamidey-noury.fr
www.lamidey-noury.fr

French manufacturer of electrosurgical unit since 1947